concept

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Digitálna dojčenská váha

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Digital baby scale


Digitale Babywaage

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Bilancia digitale per neonati

Băscula digital para bebés

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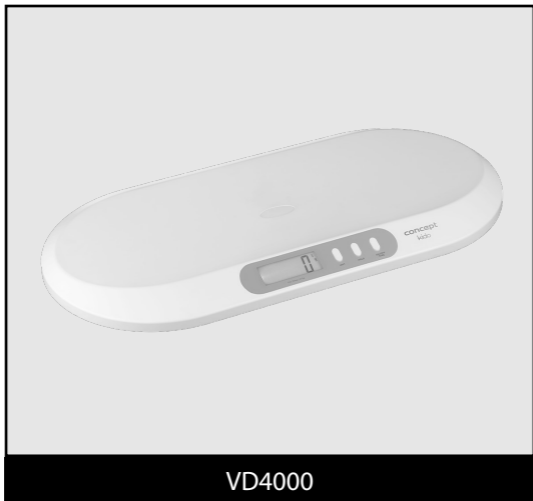VD4000

CZ

SK

PL

HU

LV

EN

DE

FR




IT

ES

RO

POPIS VÝROBKU
POPIS VÝROBKU
OPIS PRODUKTU
A TERMÉK LEÍRÁSA
IERĪCES APRAKSTS
PRODUCT DESCRIPTION

PRODUKTBESCHREIBUNG
DESCRIPTION DU PRODUIT
DESCRIZIONE DEL PRODOTTO
DESCRIPCIÓN DEL PRODUCTO
DESCRIEREA PRODUSULUI

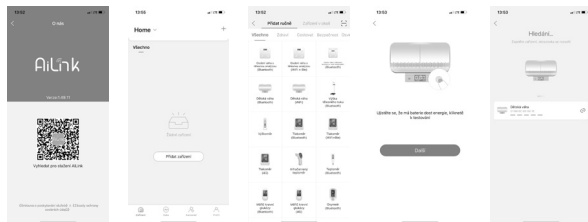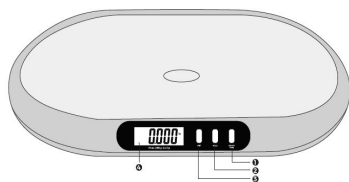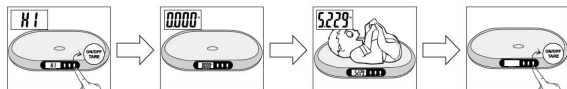Obr. 1
Fig. 1

Obr. 2
Fig. 2

Obr. 3
Fig. 3

Obr. 4
Fig. 4

Obr. 5
Fig. 5

PODĚKOVÁNÍ

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ACKNOWLEDGMENT

Thank you for purchasing a Concept product. We wish you much joy with your new appliance every day you use it.

Please read the instruction manual carefully before initial use. To refer to this manual any time you need to, we recommend you to keep it in a safe place. And please pass it to any future owner of the appliance.

TECHNICAL SPECIFICATIONS

Voltage Battery 3x AAA

IMPORTANT SAFETY PRECAUTIONS

- Do not use the appliance in any other way than as described in this manual.
- Remove all the covering and marketing materials from the appliance before the first use.
- Make sure the mains voltage corresponds to the values on the rating plate of the appliance.
- The scale should be placed on an even, dry and hard surface.
- Do not overload the scale, maximum load is 20 kg.
- Avoid shocks or impacts on the scale.
- Do not immerse the scale in water.
- When stored, no items are allowed on the scale.
- Optimal operation and storage temperature is 5°C to 30°C.
- Do not leave the appliance unattended when it is switched on or connected to the mains.
- Do not allow children or unauthorized people to handle the appliance, use it out of their reach.
- Persons with reduced mobility, reduced sensory perception, mental disadvantage or persons not familiar with the operation must use the appliance only under the supervision of responsible person, who is familiar with the instructions.
- Pay extra attention if the appliance is used near children.
- Do not use the appliance as a toy.
- Do not let the power cord touch hot surfaces.
- Do not allow the appliance to hang on the power cord.
- Do not use the appliance in an environment where there are explosive gases and flammable substances (solvents, varnishes, adhesives, etc.).
- Do not use coarse and chemically aggressive substances to clean the appliance.

- Do not use the appliance if it is not working properly, has been dropped, damaged or immersed in liquid. Have it checked and repaired by an authorized service centre.
- Do not use the appliance with a damaged cord or plug and have it repaired immediately by an authorized service centre.
- Do not wrap the power cord around the appliance body.
- Do not use the appliance outdoors.
- The appliance is only suitable for home use and is not intended for commercial use.
- **Do not immerse the power cord, plug or appliance in water or any other liquid.**
- Do not use any accessories other than those recommended by the manufacturer.
- Do not repair the appliance yourself. Contact an authorized service centre.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Cleaning and maintenance performed by the user must not be carried out by children unless they are 8 years old and under supervision. Children aged under 8 years of age must be kept away from the appliance and its power cord. Do not let children play with the appliance.

Damages from failure to follow the manufacturer's instructions are not covered by the warranty.

PRODUCT DESCRIPTION

- 1 Switching On/Off /Zero setting
- 2 Holding the scale
- 3 Units
- 4 Display

Control knob

Error	Overloading
Lo	Low battery

OPERATING INSTRUCTIONS

Buttons

1. Press ON/OFF/TARE
 - Press this button to switch on/off the device.
 - If you use towel below baby's body, put it on the platform and press the button to set the display to "0".
2. Button HOLD
 - Press this button to hold baby's weight.
3. Button UNIT
 - Press this button to set measuring units: gram (g), kilogram (kg), pound (lb) or stone (st: lb).
4. Display (LCD)

Weighting

1. Place the scale on a solid and even surface.
2. Press "ON/OFF", "Hi" will be displayed.
3. You can start weighting, when "0,00 kg" has been displayed.
4. Put the baby or other things on the scale, the weight will be shown on the LCD display
5. Press and hold "ON/OFF" to switch off the scale (or the scale will be automatically switched off after 30 seconds of inactivity).

Installation and connection of the app

1. Download AiLink app from Appstore or Google Play. (Fig.1)
2. Switch on Bluetooth on your mobile telephone and open the app..
3. Click on Add Device and select Baby Scale (Bluetooth), (Fig. 2 and 3)
4. Switch on the scale and click to Next (Fig. 4)
5. When the scale is displayed, confirm your selection (Fig.5)
6. Now you can start using the scale

CLEANING AND MAINTENANCE

CAUTION!

To clean the surface of the appliance, use only a damp cloth, no detergents or hard objects, as they can damage the surface of the appliance! **Never clean the appliance under running water, do not rinse it or immerse it in water!**

SERVICE

Extensive maintenance or repair, which requires access to internal parts of the product, has to be performed by a qualified service centre.

ENVIRONMENTAL CONCERNS

- Packaging materials and old appliances should be recycled.
- Packaging materials may be disposed of as sorted waste.
- Dispose the plastic bags made of polyethylene (PE) of as sorted waste.

Recycling of the appliance at the end of its lifespan

This appliance is labelled in accordance with European Directive 2012/19/EU concerning used electrical and electronic appliances (Waste Electrical and Electronic Equipment - WEEE). The symbol on the product or on its packaging indicates that this product may not be treated as household waste. Instead, it should be taken to the appropriate collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. The disposal of such material has to be done according to the recycling regulations. For more detailed information about recycling of this appliance see the authorized local office or service for household waste disposal or the shop, where you purchased the appliance.

The product meets all the necessary requirements of the EU directives applicable to the specific product.

Changes in text, design, and technical specifications are subject to change without notice and we reserve the right to change them.

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